

PRDM1/BLIMP1 | Validation File

TARGET Human PRDM1/BLIMP1 alpha and beta isoforms

CLONE NAME ROS195G

DESCRIPTION Mouse monoclonal

ANTIGEN USED GST-PRDM1 (full length protein)

ISOTYPE IgG1

SPECIES REACTIVITY Human and mouse

LOCALIZATION Nuclear

POSITIVE CONTROL Tonsil

STORAGE BUFFER Tissue culture supernatant: 0.02% sodium azide

STORAGE Aliquot and store at 4C. Do not freeze



Recommended



Inconclusive



Not Recommended




Not Tested

APPLICATIONS

● | ICC | Immunocytochemistry

ROS195G is able to detect human PRDM1/Blimp1 protein in immunocytochemistry

To confirm that ROS195G mAb recognizes human PRDM1/Blimp1 protein, immunocytochemistry on frozen cytospin preparations of Flag-tagged human PRDM1 expressed in HEK293T cell line was performed. HEK-CAS9-FLAG transfected cells were used as negative control.

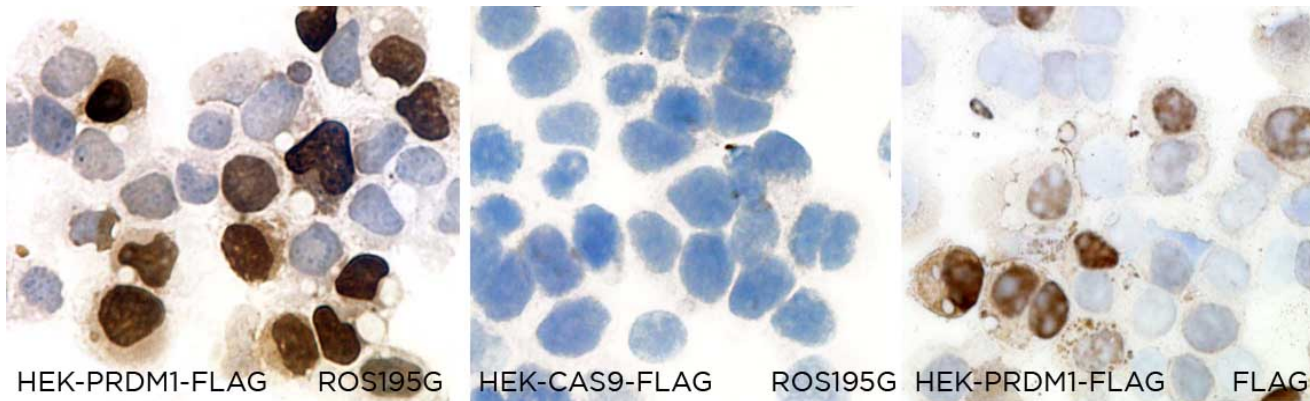


● | WB | **Western Blotting**

ROS195G mAb is able to detect human PRDM1 protein by WB.

DILUTION neat supernatant

LANES

Lane 1 L-540	(200ug) (+)
Lane 2 L-428	(200ug) (+)
Lane 3 HDLM-2	(200ug) (+)
Lane 4 KMH2	(200ug) (-)
Lane 5 L1236	(200ug) (-)
Lane 6 Human tonsil	(200ug) (+)
Lane 7 RPMI-8226	(200ug) (+)
Lane 8 KARPAS 422	(200ug) (+)
Lane 9 PBL (EBV)	(200ug) (+)

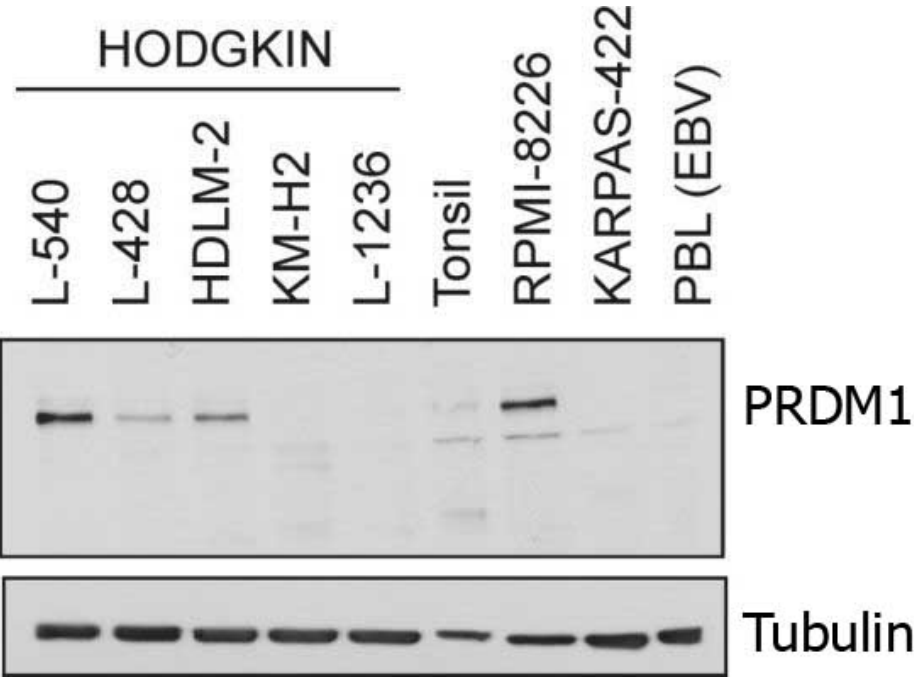


● | IP | Immunoprecipitation

Immunoprecipitation of protein extracts from normal human tonsil with ROS195G antibody and a specific polyclonal antibody (kindly provided by Dr. Mark M. Davis, USA) followed by WB with ROS195G antibody.

DILUTION Mouse antibody ROS195G was used as undiluted supernatant (40ul/lane) and anti-PRDM1 polyclonal was used 1:150.

RESULTS:

Immunoprecipitation of protein extracts from normal tonsil with the anti-PRDM1 ROS195G antibody and with one anti-PRDM1 polyclonal followed by western blotting with the ROS195G antibody (neat). In both cases it is possible to observe a specific band of 95kDa confirming antibody specificity.

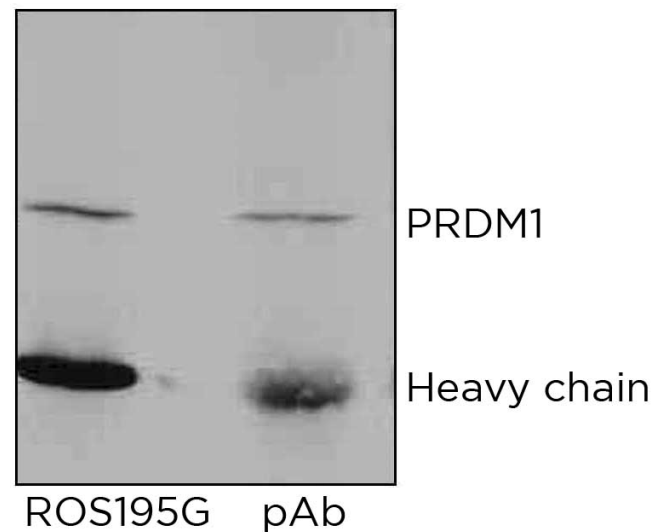


● | IHC-F | **Immunohistochemistry (frozen)**

ROS195G mAb can be used to detect PRDM1/Blimp1 protein in human frozen tissues.

TISSUE SAMPLE Human tonsil

DILUTION neat supernatant

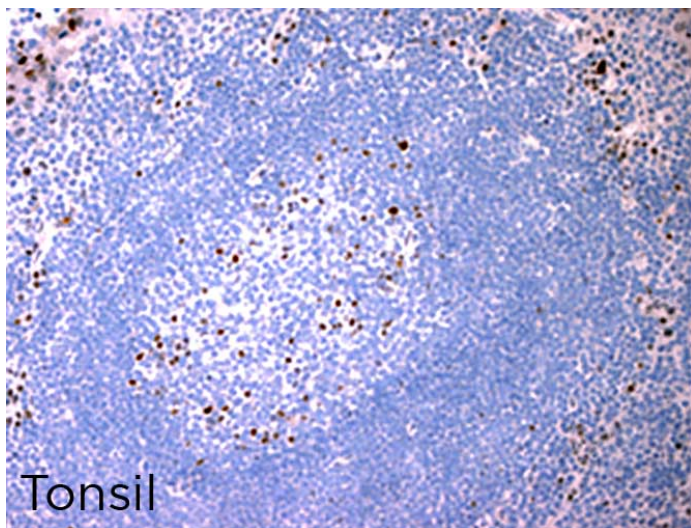


● | M-IHC-P | **Immunohistochemistry (paraffin)**

ROS195G mAb can be used to detect Blimp1 protein in mouse human paraffin tissues

TISSUE SAMPLE mouse lymph node, intestine, brain and spleen

DILUTION Neat supernatant

DETECTION SYSTEM Discovery XT (Ventana) CC1 OmniMap

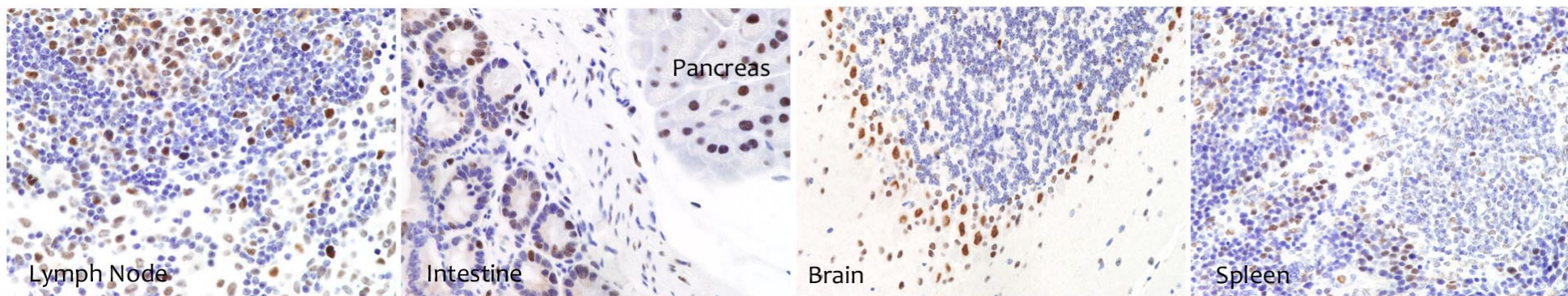

● | IHC-P | **Immunohistochemistry (paraffin)**


ROS195G mAb can be used to detect PRDM1 protein in human paraffin tissues

TISSUE SAMPLE Human tonsil and bone marrow myeloma

DILUTION 1:10

ANT. RETRIEVAL 20 minutes ER2 (Tris-EDTA)

DETECTION SYSTEM Novolink kit (BondMax Leica)

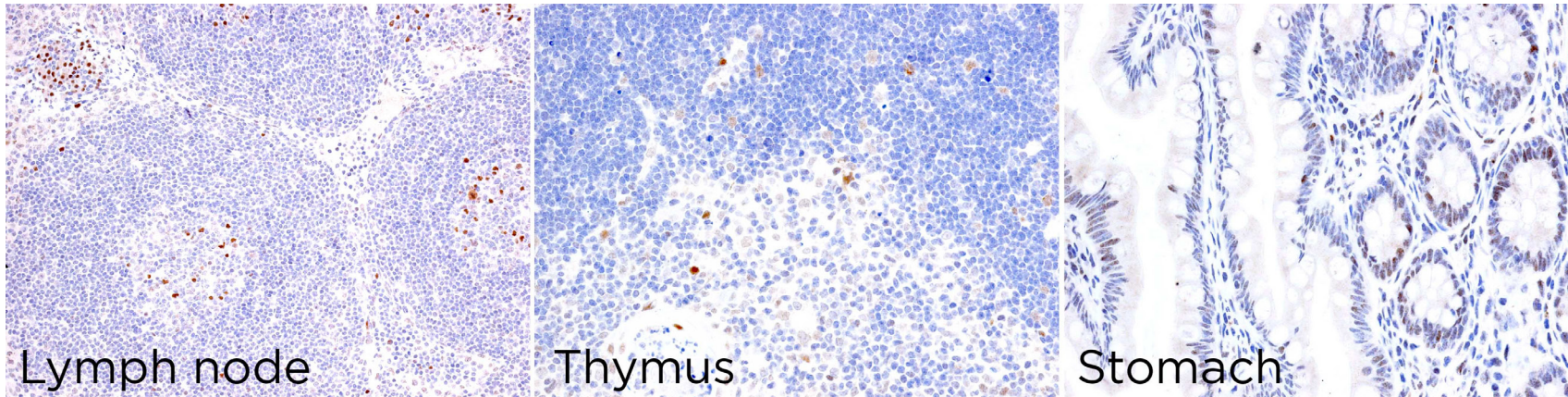

● | FC | **Flow Cytometry**


ROS195G mAb can be used to detect PRDM1 protein by Flow Cytometry

TISSUE SAMPLE human lymphoid cell lines (MOLT4 and OPM2)

DILUTION 50ul ROS195G supernatant/one million cells/tube

RED: Negative cell line MOLT4 (T cell leukemia)

BLUE: Positive cell line OPM2 (multiple myeloma)


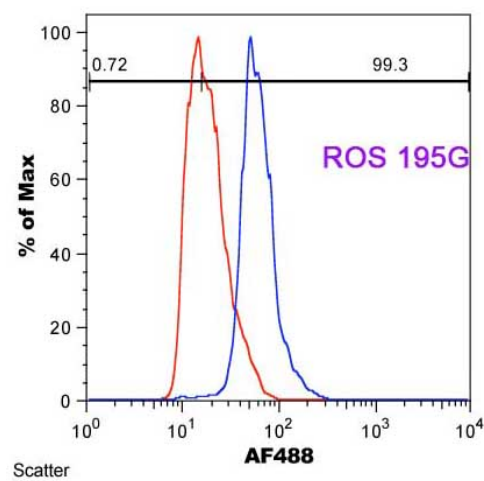


● | IF | **Immunofluorescence (paraffin)** Not Recommended

SOLD BY: Abcam and Active Motif

Additional information

ROS195G monoclonal antibody recognizes the site between amino acids 176 and 307 and is recognizing the two PRDM1/Blimp-1 isoforms. Work done by Dr K. Wright (Lee Moffitt Cancer Center).

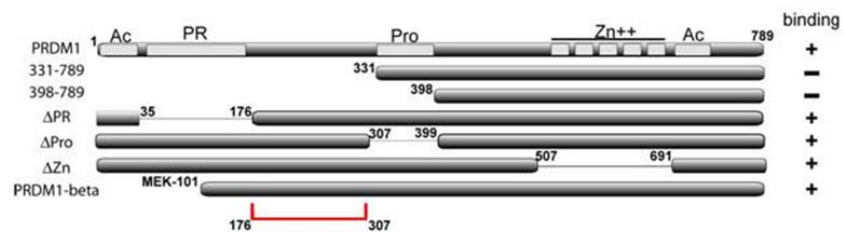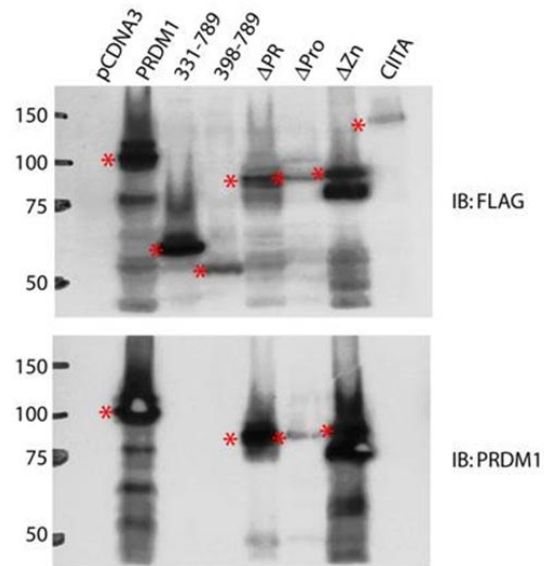

REFERENCES

Hsiao-Wei Tsao, Tzong-Shyuan Tai, William Tseng, Hui-Hsin Chang, Roland Grenningloh, Shi-Chuen Miaw, I-Cheng Ho. Ets-1 facilitates nuclear entry of NFAT proteins and their recruitment to the IL-2 promoter. Proceedings of the National Academy of Sciences of the United States of America. 24 September 2013.

Marion Travert, Yenlin Huang, Laurence de Leval, Nadine Martin-Garcia, Marie-Helene Delfau-Larue, Françoise Berger, Jacques Bosq, Josette Brière, Jean Soulier, Elizabeth Macintyre, Teresa Marafioti, Aurélien De Reyniès, Philippe Gaulard. Molecular features of hepatosplenic T-cell lymphoma unravels potential novel therapeutic targets. Blood. 14 June 2012.

Martin-Arruti M, Vaquero M, Díaz de Otazu R, Zabalza I, Ballesteros J, Roncador G, García-Orad A. Bcl-2 and BLIMP-1 expression predict worse prognosis in gastric diffuse large B cell lymphoma (DLCL) while other markers for nodal DLBCL are not useful. Histopathology. 2012 Apr;60(5):785-92

Liu YY, Leboeuf C, Shi JY, Li JM, Wang L, Shen Y, Garcia JF, Shen ZX, Chen Z, Janin A, Chen SJ, Zhao WL. Rituximab plus CHOP (R-CHOP) overcomes PRDM1-associated resistance to chemotherapy in patients with diffuse large B-cell lymphoma. Blood 2007. Jul1; 110(1): 339-44.

Zhao WL, Liu YY, Zhang QL, Wang L, Leboeuf C, Zhang YW, Ma J, Garcia JF, Song YP, Li JM, Shen ZX, Chen Z, Janin A, Chen SJ. PRDM1 is involved in chemoresistance of T-cell lymphoma and down-regulated by the proteasome inhibitor. Blood 2008. Apr 1;111(7): 3867-71.

Nie K, Gomez M, Landgraf P, Garcia JF, Liu Y, Tan LH, Chadburn A, Tuschl T, Knowles DM, Tam W. MicroRNA-mediated down-regulation of PRDM1/Blimp-1 in Hodgkin/Reed-Sternberg cells: a potential pathogenetic lesion in Hodgkin lymphomas. Am J Pathol. 2008. Jul; 173(1): 242-52.

Garcia JF, Roncador G, Garcia JF, Sanz AI, Maestre L, Lucas E, Montes-Moreno S, Fernandez Victoria R, Martinez-Torrecedrara JL, Marafioti T, Mason DY, Piris MA. PRDM1/BLIMP-1 expression in multiple B and T-cell lymphoma. Haematologica. 2006. Apr; 91:467-74.

Nam-Cha SH, Montes-Moreno S, Salcedo MT, Sanjuan J, Garcia JF and Piris MA. Lymphocyte-rich classical Hodgkin's lymphoma: distinctive tumor and microenvironment markers. Modern Pathology (2009) 22, 1006-1015.